

MicroFlow III Ductless Workstation

The MicroFlow III Workstation ideal for: A class 1 ductless hood workstation Equipped with activated carbon filtration, for fumes, odors, non-hazardous chemical vapors and pre-filter. Designed to reach into with hands and forearms only.


- Easily moved and Completely self contained with integral recessed work surface to contain spills.
- Clear hood surround with safety viewing sash allows visibility from 3 sides, and provides both safety and access to inside work area and protect vision. Sash can be conformed for use with a microscope.
- Variable speed fan control. High speed for 100f/m. Air flow thru viewing sash closed.
- Operates on 115V AC, or 230V International, conforms to UL, CSA, & CE requirements.
- Electrical cord port exits both sides.
- Easy change out of used filter. Recommend filter change every 60-90 days. Filter life varies with usage.
- Filter is removed and replaced inside the fume chamber for user safety.
- Dimensions 24" wide X 20.75" high X 24" deep. Weight allows for easy access with hands and forearms up to elbows.



Rear view of MicroFlow III

Description	Cat. No.
1. MicroFlow III Workstation 115V	74010
2. Replacement carbon filter	52006
3. Replacement pre-filter	74011
4. Mobile Table 30" wide X 24" Deep X 30" high	50303

Easy Filter Change

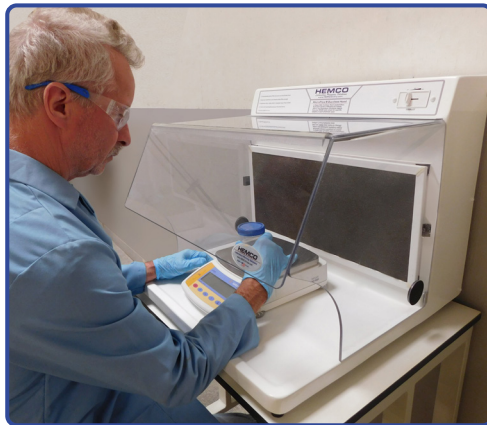

MicroFlow III



Remove plastic wrap from carbon filter & snap into place using filter clips



Install particle prefilter and attach using velcro

