

UNIMAX Floor-Mounted Walk-In Fume Hoods

The Unimax Walk-in Fume Hood encloses your workspace to safely contain and vent vapors or contaminants. Available in standard or explosion proof models in 72", 96", 108" or 144" widths



UNIMAX Floor-Mounted Hood #50601 shown with optional table, apparatus grid, and service fixtures.

Built for Durability and Convenient Use

Manufactured of chemical resistant composite materials bonded to 2" urethane foam. Horizontal sliding doors are constructed of insulated tempered glass. The Unimax is available in "double depth" or "walk-through" designs to provide convenient front and rear access. Unit can be mounted to existing floor or to other special floor covering as required.

Air By-pass Design Controls Air Flow


The automatic Air-By-Pass system controls air velocities. A chemical resistant adjustable and removable slotted baffle provides uniform air flow throughout the fume chamber.

Custom Designed to Customer's Needs

Equipped with either recessed vapor light or explosion proof incandescent light. Chemical resistant exhaust blowers for normal or explosion proof applications are available. Can be shipped completely assembled or "knocked down" for on-site assembly.


Distillation Apparatus Grid

Offered in aluminum, stainless steel, and non-metallic composite plastic. Includes rod and all connectors. The 76" hood requires one (1) 6' assembly, the 100" hood requires two (2) 4' assemblies, and the 148" hood requires three (3) 4' assemblies.

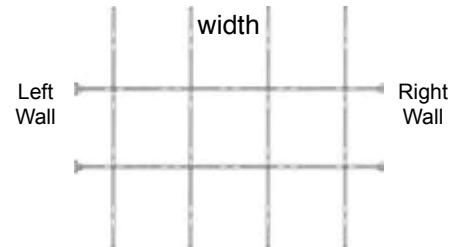

	76" Width	100" Width	112" Width	148" Width
	Qty Cat No	Qty Cat No	Qty Cat No	Qty Cat No
1. Aluminum Rod Assy.	(1) 50261	(2) 50251	(2) 50251	(3) 50251
2. Stainless Steel Assy.	(1) 50262	(2) 50252	(2) 50252	(3) 50252
3. Plastic / FG Rod Assy.	(1) 50263	(2) 50253	(2) 50253	(3) 50253

UniMax Floor Mount Fume Hoods

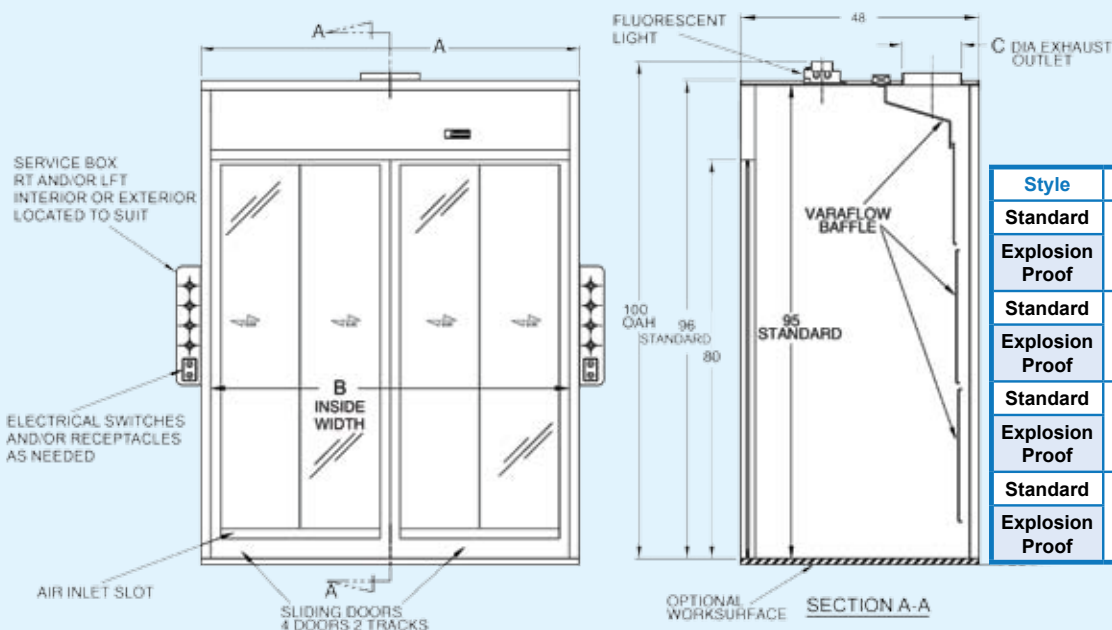

Style	A	B	C	Cat. No.
Standard				50601
Explosion Proof	76"	72"	12"	50631
Standard				50602
Explosion Proof	100"	96"	12"	50632
Standard				50603
Explosion Proof	112"	108"	12" (2pl)	50633
Standard				50604
Explosion Proof	148"	144"	12 (2pl)	50634

Also available with access doors on front and back. For plumbing, sinks and electrical fixtures, see pages 27-29. Consult HEMCO for custom UNIMAX applications