

Student Demonstration Hood (Ductless)

The Ductless Demo Hood effectively removes fumes and odors while providing a clear, unobstructed view of the interior from all sides. The demonstration hood is designed for classroom and group demonstrations.




- A Class DH-1 ductless hood workstation equipped with activated carbon filtration to control non-toxic chemicals, fumes and nuisance odors.
- Dimensions 36" or 48" wide X 24" deep X 30.5" high.
- Viewing height is 24", through hinged clear acrylic viewing sash, shatterproof side and rear panels.
- Front control panel includes light switch, and exhaust switch, hood comes equipped with a vapor proof T-5 fluorescent light fixture, 158 CFM blower, and a 8' power cord.

Specifications:

Ductless Demonstration Hood: 115 V AC, 60 Hz, shipping weight 250 lbs	Cat. No. 53015
Demo Hood equipped for International Electrical 230 V AC, 50 Hz	Cat. No. 53015
Height Adjustable Mobile Table: Height adjustable from 25" to 37", Top measures 24" X 36"	Cat. No. 50318
Non-Adjustable Table Fixed height 36"	Cat. No. 50315
Replacement Filter: Qty 2 pack	Cat. No. 52005



Easy Filter Change

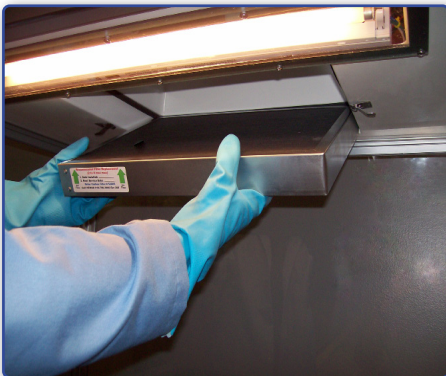

Release the filter clips and remove old filter



Install new filter and re-latch filter clips



Demo hood is ready for use

Visit www.HEMCOcorp.com

