


Model 14000 Table Top Hood

Portable Table top Canopy hoods provide convenient light duty fume removal and are ideal as multiple stations for student or industrial use.

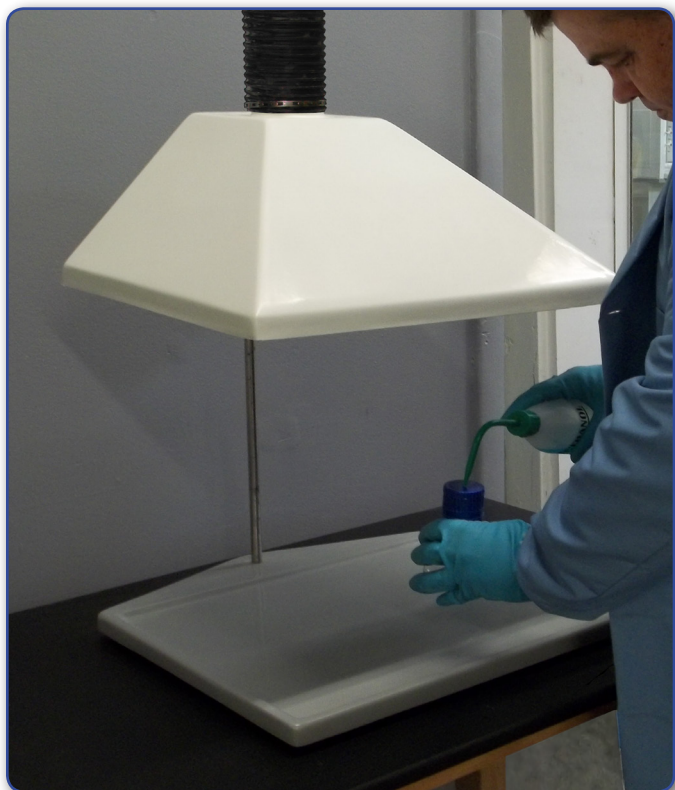


- Constructed of chemical resistant advanced composite materials, this table top hood can be adjusted in height up to 18" to accommodate your specific instrument or procedure. "Easy-Touch" operation allows for precise positioning.

- Measuring 18" wide x 18" deep and 25" tall
- Lightweight, it can easily be moved as needs change.
- Recessed work surface to contain spillage.
- Hood can easily vent out a wall, roof, or into a existing exhaust system.
- Unit features a 3" dia. outlet collar for duct connection.
- Shipped completely assembled, ready for operation.

Accessories:

Flexible Ducting <ul style="list-style-type: none">• 3" flexible PVC material includes(2) stainless steel clamps, easy to install, ships in 12' lengths	Cat No. 80047
Exhaust Blower with mounting bracket <ul style="list-style-type: none">• Mounts to wall or table. Epoxy coated steel blower housing and blower wheel, 135 CFM @ .3"sp,3" dia inlet and outlet, includes mounting bracket	Cat No. 51050
Vent Outlet <ul style="list-style-type: none">• For venting through the wall or window, 3" dia.,constructed of PVC pipe and clear acrylic mounting plate.	Cat No. 80073

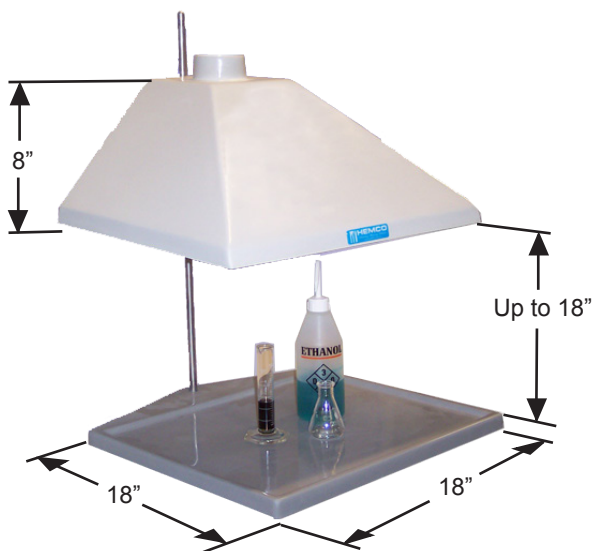

Typical Applications

- Student Workstations
 - Dental Labs
- Microscope Station
- Microchemistry
- Sample Holding
 - Pipetting
 - Microbiology

Visit www.HEMCOcorp.com