

Model 24000 Vented Table Top Hood

A vented hood workstation. HEMCO's 24000 Table Top features a molded seamless construction, including rear and side walls to minimize cross drafts in the work area.

Dimensions are 24" wide X 15" deep X 24" high


- Constructed of chemical resistant, light weight advanced composites, it can be easily moved as procedures or workflow change.
- Molded chemical resistant work surface is recessed to contain spillage, and a three inch diameter outlet collar is provided for duct connection.
- Fumes are vented through the integral fume side and rear walls and out the top. Cross drafts are reduced, thus providing a turbulence free work area.
- The units base consists of an integral recessed work surface to contain spillage.
- Typical uses include, Histology, Microprocessor, Venting for Hot Plates, Microscope Stations, Student workstations, Sample weighing stations, and handling Pharmaceuticals.



| Model 24000 Vented Table Top Hood | Cat. No. 24000 |
|--|----------------|
| Accessories | Cat. No. |
| Flexible Ducting - 3" Diameter Flex ducting with stainless steel clamps, in 12' lengths | 80047 |
| Exhaust Blower - mounts to wall or table. Epoxy coated steel housing and wheel, 135 CFM @ .3"sp, 3" dia inlet and outlet, includes mounting bracket | 51050 |
| Vent Outlet - for through the wall or window venting. 3" dia., constructed of PVC Duct with acrylic mounting plate | 80073 |
| Vapor Proof Light - factory installed with a 6" power cord and switch. 115V U.L listed | 50035 |
| Clear Hinged Safety Shield - provides extra protection for user, size 6" X 6" X 11" | 24121 |



1. 24000 with Safety Shield in the Up position

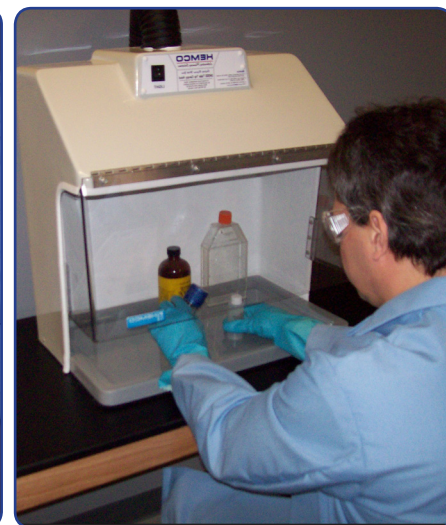

2. 24000 with Safety Shield in the Down position

Visit www.HEMCOcorp.com

HEMCO Corporation
711 S. Powell Road
Independence, MO 64056

HEMCO

Phone (816) 796-2900
Fax (816) 796-3333
24000-08.11