


NEW UniLine ACID Storage Cabinet



For storage of ACID and corrosive chemicals, available in 12", 18", 24", 30", 36", 42" & 48" widths. Standard 35" high and 22" deep cabinet dimension. (as shown) U.L. 1805 Fume Hoods & Cabinets

- The molded one piece fiberglass liner inserts directly in the cabinet and is sealed on all edges with 1 1/4" spill containment lip, across bottom front edge. Rear Access Panel
- The Front Access doors with air inlet vents are lined and edges are sealed. No metal is exposed to corrosive vapors.
- Removable Shelf for small container storage

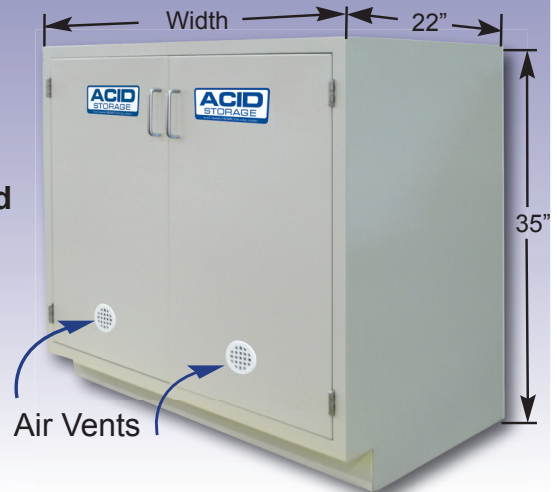

Acid Cabinet Size							
Width	48"	42"	36"	30"	24"	18"	12"
Catalog No.	15480	15035	15360	15300	15240	15180	15005
All Acid Cabinets are 22" deep X 35 1/4 high							

Venting Your Acid Cabinet (select your choice 1.2.3...)


2" Acid Cabinet Vent Through Cabinet Back



① Acid Cabinet Vent Through Fume Hood Work Surface

- (1) 90 Degree Elbow
- (1) Vent outlet Exhaust Vent Kit


Cat. No. 81000

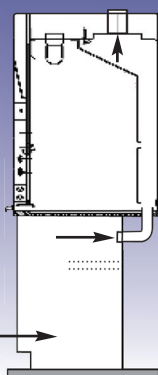

② Acid Cabinet Vent Through Fume Hood Double Side Wall

- (1) 90 Degree Elbow
- (1) Vent outlet Exhaust Vent Kit


Cat. No. 81001

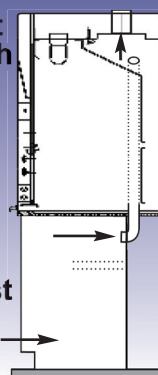

③ Acid Cabinet Vent Thru Building Wall with Vent Blower

- (1) 90 Degree Elbow Exhaust
- (1) 10' Ridge & Flex
- (1) Exhaust Blower Vent Outlet

Cat. No. 81002

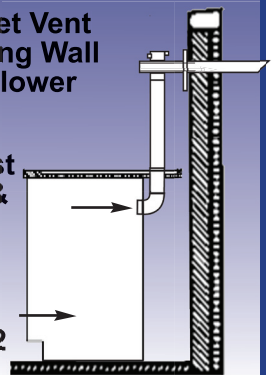

Scientific Equipment Furniture Association



Meets NFPA-45 for Product used in Laboratories



Greener Product Certified



Underwriters Laboratories UL 1805 Classified



HEMCO™

Visit www.HEMCOcorp.com