

NEW UniLine ACID Storage Cabinet



For storage of ACID and corrosive chemicals, available in 18", 24", 30", 36", & 48" widths. Standard 35" high and 22" deep cabinet dimension. (as shown)

- The molded one piece fiberglass liner inserts directly in the cabinet and is sealed on all edges with 1 1/4" spill containment lip, across bottom front edge.
- The Front Access doors with air inlet vents are lined and edges are sealed. No metal is exposed to corrosive vapors.
- Shelf for small container storage



Acid Cabinet Size					
Width	18"	24"	30"	36"	48"
Catalog No.	15180	15240	15300	15360	15480
All Acid Cabinets are 22" deep X 35 1/4 high					


Venting Your Acid Cabinet (select your choice 1.2.3...)

Acid Cabinet Vent Through Cabinet Back




① Acid Cabinet Vent Through Fume Hood Work Surface

- (1) 90 Degree Elbow
 - (1) Vent outlet Exhaust Vent Kit
- Cat. No. 81000

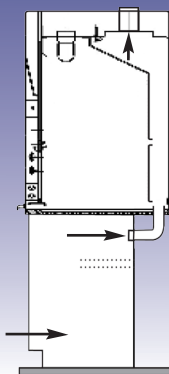


② Acid Cabinet Vent Through Fume Hood Double Side Wall

- (1) 90 Degree Elbow
 - (1) Vent outlet Exhaust Vent Kit
- Cat. No. 81001

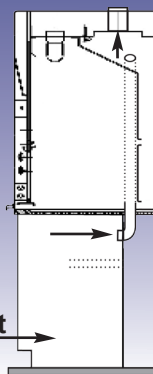

③ Acid Cabinet Vent Through Building Wall

- (1) 90 Degree Elbow Exhaust
 - (1) 10' Ridge & Flex
 - (1) Exhaust Blower
 - (1) Vent Outlet
- Cat. No. 81002

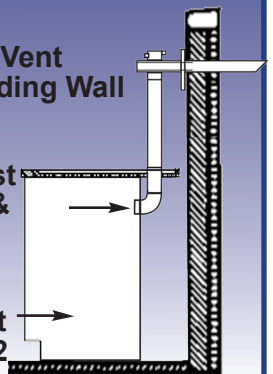

Meets NFPA-45 for Product used in Laboratories



Greener Product Certified



Scientific Equipment Furniture Association



HEMCOTM
Laboratory Planning Solutions

Visit www.HEMCOcorp.com